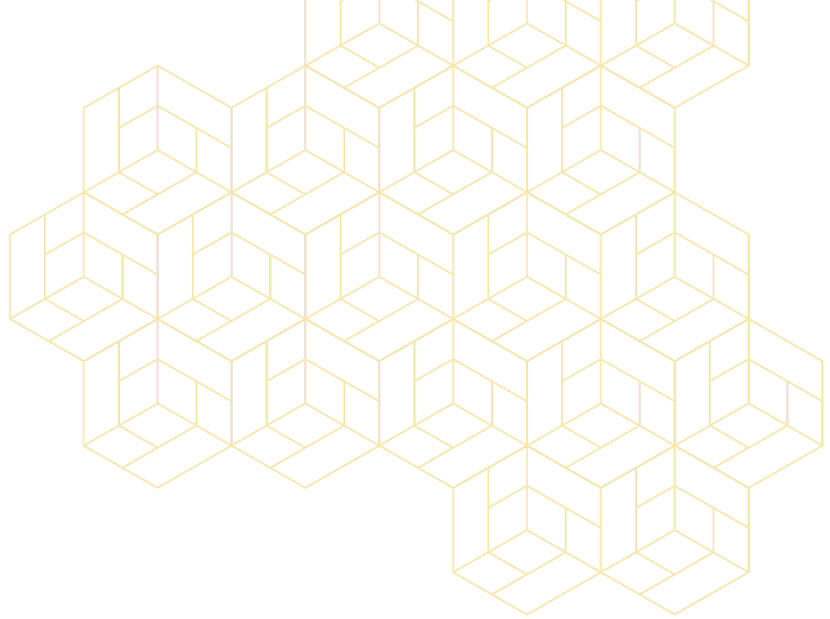




[Test Update 967](#)

Posted Date 08/12/2026

Effective Date 08/25/2026

Test Name [Hantavirus Antibodies, IgG & IgM](#)

Update Type [New Tests](#)

CPT Code 86790 x2

NEW TEST

Hantavirus Antibodies, IgG & IgM

Order Code: pending

Fee Code: AA491 (CPT 86790 x2)

Reference Laboratory: Mayo HANGM

Mayo Clinic Laboratories will begin offering Hantavirus Antibodies, IgG & IgM, effective August 25, 2026. This test will replace Mayo test ID FHVGM (order code FHVGM) referred to Quest Diagnostics.

Collection Instructions: Collect specimen in an SST (preferred) or red top tube. Centrifuge, aliquot 0.6 mL (minimum 0.5 mL) of serum into a plastic vial, and refrigerate.

Reference Range: IgM: Negative, IgG: Negative.

- False-negative results may occur in patients that are tested either too soon following infection or are significantly immunocompromised.
- False-positive results may occur in patients with recent cytomegalovirus infection.
- This assay detects IgM and IgG antibodies against antigens from the Sin Nombre and Andes Hantavirus strains; however, it will not differentiate between the two strains. It is unknown if this assay will detect antibodies to other Hantavirus strains.
- Results must be interpreted alongside epidemiological exposure history and clinical presentation, and diagnosis should not be made solely on serologic results from this assay.

[View PDF](#)