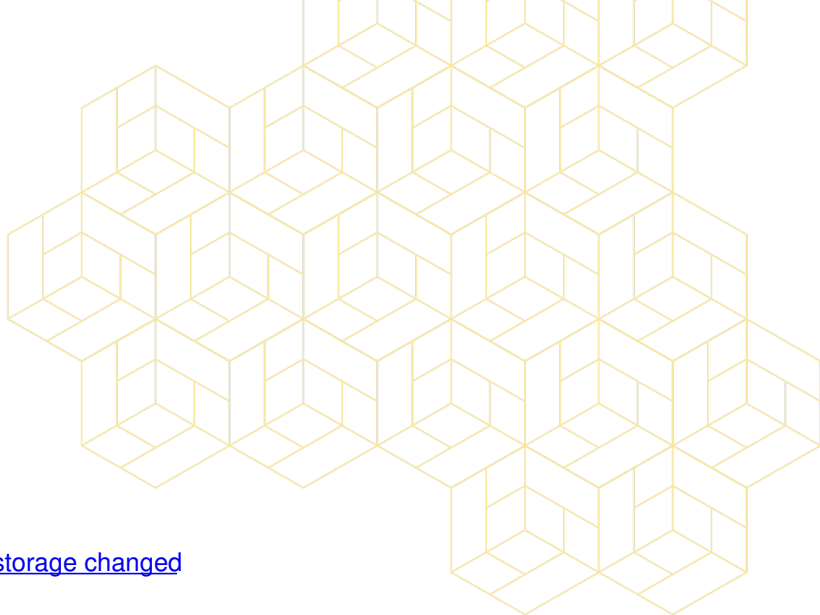




Test Update 968

Posted Date 08/19/2026

Effective Date 08/07/2026

Test Name [Zinc Transporter 8 \(ZnT8\) Antibody](#)

Update Type [Specimen collection or handling and/or storage changed](#)

SPECIMEN HANDLING CHANGE

Zinc Transporter 8 (ZnT8) Antibody

Order Code: FZNT8

Fee Code: AA634 (CPT 86341)

Reference Laboratory: Mayo FZNT8 (Esoterix Endocrinology)

Please note the follow change to Zinc Transporter 8 (ZnT8) specimen collection and handling instructions effective August 8, 2026:

Collection Instructions: Collect specimen in a red top (preferred) or SST tube. Centrifuge, aliquot 0.5 mL of serum into a plastic vial within 45 minutes of collection, and freeze.

[View PDF](#)